

# THE EMERGENCY BANDAGE

## Trauma/Wound Compression Dressing

**The first major innovation of its kind in over 50 years!**

- NSN #6510-01-460-0849 -- FCP-01, 4" width
- NSN #6510-01-492-2275 -- FCP-02, 6" width
- NSN #6510-01-515-7528 -- FCP-05, 6" + a mobile pad
- NSN #6510-01-532-6656 -- FCP-09, 8" abdominal




**Performance Systems** presents an innovative, combat proven first-aid device for the staunching of blood flow from traumatic hemorrhagic wounds in pre-hospital emergency situations.

**This all-in-one device** consolidates multiple first-aid devices such as a primary dressing, pressure applicator, secondary dressing, and a foolproof closure apparatus to secure the bandage in place. This internationally patented and FDA approved bandage is the ideal solution for emergency treatment when every second counts.

**The Emergency Bandage** is designed for quick and easy application by professional and non-professional caregivers to provide effective, multi-functional treatment.

**The Emergency Bandage's** sterile, non-adherent pad applies pressure to any site, can be easily wrapped and secured, and has an additional application, similar to a tourniquet, to further constrict blood-flow. The product is so easy to use that even an injured person can self-apply the bandage with one hand.

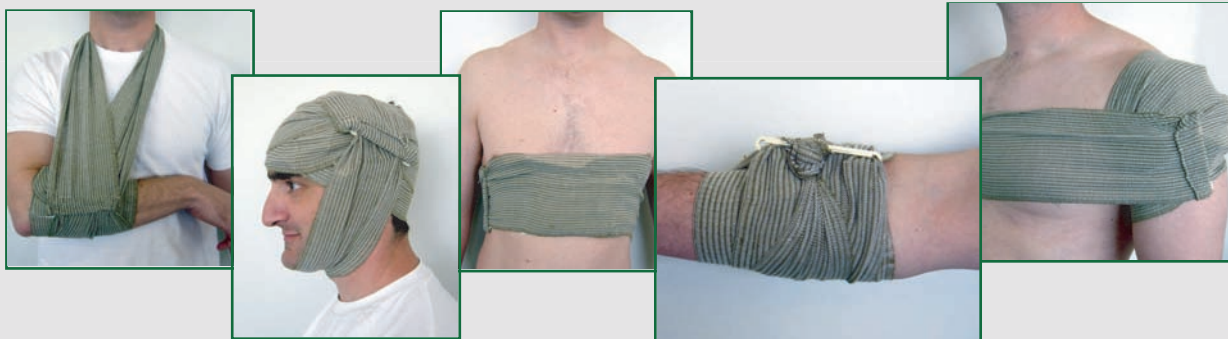

### Applications:

- Military - Medics, Soldiers
- First Responders
- Emergency Medical Services
- Law Enforcement
- Disaster Preparedness
- Homeland Security
- Commercial and Industrial Sites
- First-Aid Kits

### Advantages:

- Immediate direct pressure
- Quick and easy self-application
- Consolidation into a single unit: *primary dressing, pressure applicator, secondary dressing, closure bar*
- Significant time, space, and cost savings
- X-Ray friendly
- Versatile
- Secure, water-resistant closure bar




### Performance Systems, Medical Division

3050 Post Oak Blvd., Ste. 1710, Houston, TX, 77056  
 Toll Free 1-888-737-7978 | Tel (713) 723-6000 | Fax (713) 723-6221  
 E-mail: [info@performancesystems.com](mailto:info@performancesystems.com) | [www.ps-med.com](http://www.ps-med.com)

**NEW!**

# Emergency Bandage with Second, Mobile Pad FCP-05 Abdominal Bandage FCP-09

**FCP-09**

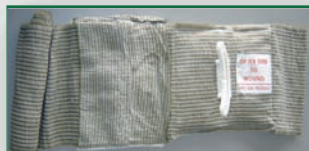

A second, mobile pad

**FCP-05**



For entry and exit wounds



For extra absorbency

## Testimonials

### **1LT Mike Nesbitt, PA-C Battalion, Physician Assistant, 1-9 FA Battlekings**

"Thank you very much for your effort to get the EB to us just before Op Iraqi Freedom.

During the war the bandage performed spectacularly. Every medic in my unit had them in his aid bag, we used it in our combat lifesaver bags and our aid station used it as the primary dressing for bleeding wounds. It performed well and was easy to teach to the non-medical personnel. It was less expensive to equip my medics with your dressing than the standard equipment, and took up less space in their aid bags. My chain of command is sold on your dressing...we will continue to use it as my unit's primary trauma dressing, replacing the standard field dressing."

### **Carl Bernstein, HM1/ 18D / SEAL TEAM 7**

"We just debriefed from Seal Team 3, your bandage was a big hit in Afghanistan. Our Seal corpsmen used them on Seals, Army and Marine wounded. All our Seal platoons will now deploy with your dressing!"

### **SFC Robert M Miller, Medical Plans, Requirements and Standardization, 75th Ranger Regiment**

"As of 10 December 2002 the Emergency Bandage is the official primary battle dressing for the 75th Ranger Regiment, United States Special Operations Command. The EB was selected as the "Best in Class" based on manufacturing, effectiveness, ease of use and price."

### **Cpl. Peter Roehrig, Pasco Sheriff's Office**

"Sept. 11, 2001 will always be a day in history no one will ever forget. Now Sept. 12, 2002 will live on as a day in history for the Pasco Sheriff's Office because on that day while executing a narcotics search warrant one of my brother deputies was shot in the head just above the right eye, penetrating and shattering his sinus cavity. Two weeks prior to this incident I attended a Tactical Medic Class, at this location I learned about your bandage and its amazing features. On the evening of Sept. 12, 2002 I placed one of the sample bandages in my SWAT Medic bag and once I assessed the Deputy's wound I by-passed any traditional bandages and immediately used the EB because of its effectiveness in controlling and securing to head wounds in a matter of seconds. The Deputy spent only four days in the hospital and was released pending further evaluation, but your bandage was a Godsend for me that night. I am now working on issuing the SWAT team these bandages as standard equipment."

### **William Madden, Medical Operations, 10th WMD CST**

"Those exposed to the bandage have been overwhelmingly amazed and interested in its ease of use and function. Especially when compared to the "traditional" field pressure dressing. I found it interesting just how quickly "non medical" personnel were able to adapt to the bandage and use it properly."

### **Bob Scavelli, PA - C**

"The emergency bandage has proven worth its weight in gold during combat maneuvers in support of Operation Iraqi Freedom. I am really impressed with the versatility of the bandage combined with the speed and ease of use. During use in treatment of patients with multiple wounds, including the unique tourniquet application, the emergency bandage has saved time and lives. The emergency bandage is well on its way to being one of the greatest assets emergency personnel can have."

## Instructions for Use

**1**



Place pad on wound and wrap the elastic bandage around limb or body part

**2**



Insert elastic bandage completely into pressure bar

**3**



Pull elastic bandage back over top sides of pressure bar forcing it down onto the pad

**4**



Wrap elastic bandage tightly over pressure bar and wrap over all edges of pad

**5**



Secure hooking ends of closure bar into elastic bandage